

Axygen[®] Long Reach Pipet Tips

Ideal for pipetting from long and narrow tubes



When pipetting, are you currently disturbed because you are unable to aspirate precious reagents/samples from long and narrow storage tubes? Does it bother you that you might contaminate the micropipettor during pipetting?

Your help is here. Axxygen[®] Long Reach Tips are elongated, clear, and compatible with vast majority of research grade pipettor brands. The tips are suitable for basic research experiments, and also able to meet the needs of high-end precision pipetting.

Features

- ▶ Slender Elongated Tip
- ▶ Graduation marks for easy visualization of liquid volume
- ▶ Detectable DNase/RNase/Human DNA free, Non pyrogenic
- ▶ Autoclavable (121 °C, 15 minutes)
- ▶ Available in various packaging sizes
- ▶ Products are manufactured in a facility that conforms to ISO14644-1 class 8 ISO requirement.
- ▶ Materials conforming with US FDA USP Class VI level
- ▶ Resin is delivered by central control feeding system during the manufacturing process, hence eliminating external contamination

Using the 10 μ l long Reach Tips, you can pipette antibody (100 μ l / tube, 1mL tube) directly from the bottom of the tube. Additionally, because the tips are longer and thinner with a smaller surface area, this reduces the volume of residual antibody that sticks on the outside of the pipette tip compared to using normal tips.

The length of 200 μ l long reach tip ensures that the micropipettor shaft and tip ejector arm is not contaminated with dyes during buffer loading or cell staining procedures.

If you work with Matrigel, we recommend the use of 200 μ l or 1mL extended tips. The tips can reach the bottom of the Matrigel bottle. The shaft and ejector arm of the pipettor will not come into contact with the mouth of the bottle and hence reducing the chances of contamination.



Figure 1. Changing the media of T25 cell culture flask with a 1 mL Long Reach Tip is simple, efficient and convenient


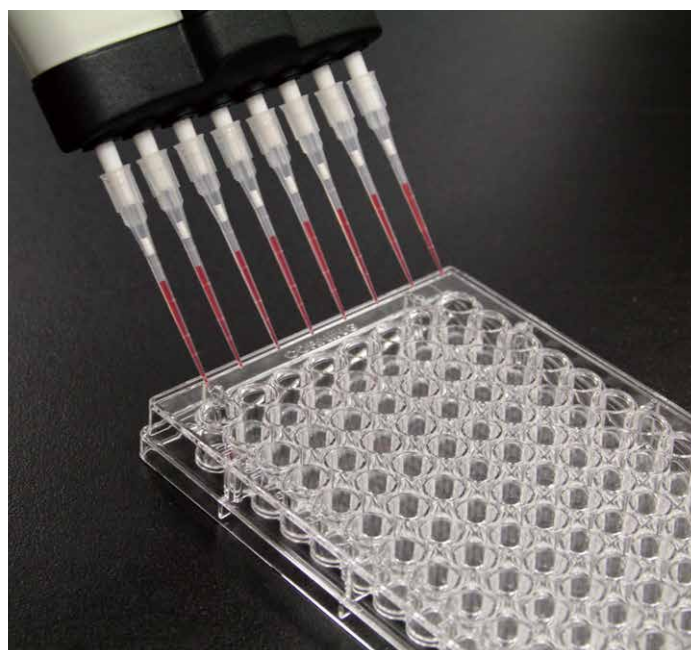


Figure 2. Micro-pipetting with 10 μ l Long Reach Tips for precision pipetting

10 µl Long Reach Pipet Tips

- ▶ Up to 45mm, longer than standard pipette tips, suitable for micro-pipetting, allows easy access to the bottom of tubes
- ▶ Graduation marks at 2µl, 5µl, and 10µl
- ▶ Volume Range: 0.1µl-10µl
- ▶ Universal fit and compatible with majority of micropipettor brands

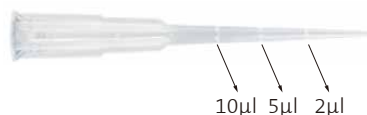


200 µl Long Reach Pipet Tips

- ▶ Up to 77.5mm, you can easily reach the bottom of the container, without damaging the pipetting pipette. 10µl, 50µl, 100µl and 200µl Minimum fluid volume 1 µl Compatible with the vast majority of pipettes
- ▶ Graduation marks at 2µl, 5µl, and 10µl
- ▶ Volume Range: 1µl-200µl
- ▶ Universal fit and compatible with majority of micropipettor brands



1000 µl Long Reach Pipet Tips

- ▶ Extended length (102mm) reaches the bottom of tubes without contaminating the pipettor shaft
- ▶ 200µl, 500µl and 1000µl mark at the mark Compatible with the vast majority of pipettes Deep hole containers are most suitable
- ▶ Universal fit and compatible with majority of micropipettor brands
- ▶ Ideal for pipetting from deep vessels

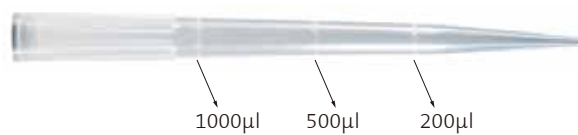

Ordering Information

CAT.NO.	Description	Packaging Format
T-10XT-C	10µl Long Reach Tips, Clear, Bulk Pack, Nonsterile	1000 Tips/Pack, 10 Packs/Case
T-10XT-C-R-S	10µl Long Reach Tips, Clear, Rack Pack, Sterile	96 Tips/Rack, 10 Racks/Pack, 5 Packs/Case
TF-10XT-C-R-S	10ul Filter Long Reach Tips, Clear, Rack Pack, Sterile	96 Tips/Rack, 10 Racks/Pack, 5 Packs/Case
T-200XT-C	200µl Long Reach Tips, Clear, Bulk Pack, Nonsterile	1000 Tips/Pack, 10 Packs/Case
T-200XT-C-R-S	200µl Long Reach Tips, Clear, Rack Pack, Sterile	96 Tips/Rack, 10 Racks/Pack, 5 Packs/Case
TF-200XT-C-R-S	200ul Filter Long Reach Tips, Clear, Rack Pack, Sterile	96 Tips/Rack, 10 Racks/Pack, 5 Packs/Case
T-1000XT-C	1000µl Long Reach Tips, Clear, Bulk Pack, Nonsterile	768 Tips/Pack, 10 Packs/Case
T-1000XT-C-R-S	1000µl Long Reach Tips, Clear, Rack Pack, Sterile	96 Tips/Rack, 8 Racks/Pack, 5 Packs/Case
TF-1000XT-C-R-S	1000ul Filter Long Reach Tips, Clear, Rack Pack, Sterile	96 Tips/Rack, 8 Racks/Pack, 5 Packs/Case



in New Zealand please contact
 Bio-Strategy Limited
 T: 0800 34 24 66 | E: sales.nz@bio-strategy.com
 www.bio-strategy.com | shop.bio-strategy.com

CORNING

Corning Incorporated
 Life Sciences
 836 North St.
 Building 300, Suite 3401
 Tewksbury, MA 01876
 t 800.492.1110
 t 978.442.2200
 f 978.442.2476

www.corning.com/lifesciences

Technical Support Hotline

Singapore: +65 6622 3609
 Australia: +61 2 8098 9173
 New Zealand: +64 9913 5888
 ScientificSupportAP@corning.com

ASIA/PACIFIC

Australia/New Zealand
 t 61 427286832

Korea
 t 82 2-796-9500
 f 82 2-796-9300

Singapore

t 65 6572-9740
 f 65 6735-2913

India

t 91 124 4604000
 f 91 124 4604099