

Newest Range Of FT/IR Instruments: The X-Series

A Comprehensive Range of FT/IR Spectrometers

(Fourier Transform Infrared Spectrometer)

The FT/IR-4X is a powerful Mid-IR FTIR spectrometer, with many features that you find in a research grade instrument, such as non-hygroscopic KRS-5 windows to prevent damage to the interferometer, a temperature controlled DLaTGS detector and a high output ceramic source for maximum sensitivity.

Permanent optical alignment is guaranteed by corner cube (retro-reflective) mirrors with auto-alignment for maximum energy.

A standalone instrument that is compatible with a wide range of sampling accessories like standard ATR, to microscopes for imaging & chemical mapping.

Please inquire for more information about alternate models in the X-Series range, configuring your instrument, or multi-unit purchasing.

Introductory Offer

\$28,600

Usually \$33,370

Spectrometer with Spectra Manager 2.5, including KnowItAll spectral database

Code: JASC7183-J006AV



[watch >](#)


Accessories: *available at additional cost*



Large Sample Compartment



Options for FTIR Microscope
and IR Imaging

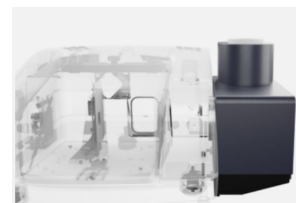

Second Detector



Prices and discounts may be subject to change and only available while stocks last. Not valid with any other special offer or discounts and cannot be applied in retrospect. All prices in local currency excluding GST and local delivery charges. Products advertised may be for research use only (RUO) or in vitro diagnostic use (IVD). Please refer to manufacturer's product insert for further information. For product use outside of that specified on the manufacturer's labelling and/or instructions for use, Bio-Strategy does not take responsibility.

Ref. Code:
JASC-FT/IR-4X

VALID UNTIL:
15 Dec 2023