

TwinGuard® ULT 729L Freezer

Only \$19,995

Was \$26,500

Code: MDF-DU702VX-PE


pricing valid for current stock only

Included **FREE**

6 x stainless steel inventory racks, (6H x 4D) with
6 pull out drawers for 2" cryoboxes
(value \$2,100.00)

The TwinGuard® -86°C ultra-low temperature freezer offers 729L of storage safety due to the dual cooling system. Two independent refrigeration systems provide the best redundancy you can have! – if one refrigeration circuit fails – the other circuit will maintain -70°C.

- Quiet in operation – only 52dB(A)
- Filter-less condenser design: eliminates the need for routine filter cleaning
- Storage for 576 x 2" cryoboxes
- Automatic vacuum relief port for quicker access to samples after door closure
- Hot gas bypass combined with fivefold door gasket assists in minimising frosting and ice build up around door seal
- Alarms for high and low temperature, door open and power failure
- Three access ports
- Intelligent ECO Mode: optimising compressor operation to match the freezer's running condition and safe electricity costs
- User friendly Colour Graphical Control Panel with digital display
- USB port for downloading logged data
- VIP Plus insulation: provides typically 24% more storage capacity compared to conventionally insulated freezers with the same footprint
- Castors for mobility and locating feet when in working position
- Keyed door locks for security
- Dry remote alarm contacts
- Suitable for operation on 220/240V 50Hz 10 amp supply
- 5-year warranty

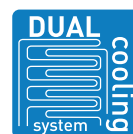

Prices and discounts may be subject to change and only available while stocks last. Not valid with any other special offer or discounts and cannot be applied in retrospect. All prices in local currency excluding GST and local delivery charges. Products advertised may be for research use only (RUO) or in vitro diagnostic use (IVD). Please refer to manufacturer's product insert for further information. For product use outside of that specified on the manufacturer's labelling and/or instructions for use, Bio-Strategy does not take responsibility.

Ref. Code:
SANY-Twin729

VALID UNTIL:
31 Mar 2024