

165L CO₂ Incubator

Only \$11,950

Was \$14,059


Code: MCO-170AICUVDL-PE

valid for current stock only

Superior culturing with the ease of reactive contamination control via high temperature sterilisation and pro-active contamination control via the inCu-saFe corrosion-resistant stainless steel copper enriched alloy interior!

- Dual heat 180°C high temperature sterilisation
- Inner glass door for viewing cultures
- Temperature: Ambient +5°C to 50°C (±0.1°C)
- CO₂: 0 to 20% (±0.15%)
- Humidity: 95%RH (±5%RH)
- Direct heat and air jacket for environment stability
- Dual wavelength single beam IR CO₂ sensor with continuous auto zero calibration
- 4 adjustable shelves
- Programmable SafeCell UV system (isolated from cell cultures) irradiates airflow and humidity pan to prevent contamination without affecting cell cultures in-vitro
- Field reversible door
- Peltier powered dew stick to draw condensation away from inner glass door
- PID control system with colour LCD display
- USB port for data transfer
- Fully rounded corners and integrated shelf supports reduce contamination and simplify cleaning
- 5-year warranty

In Stock

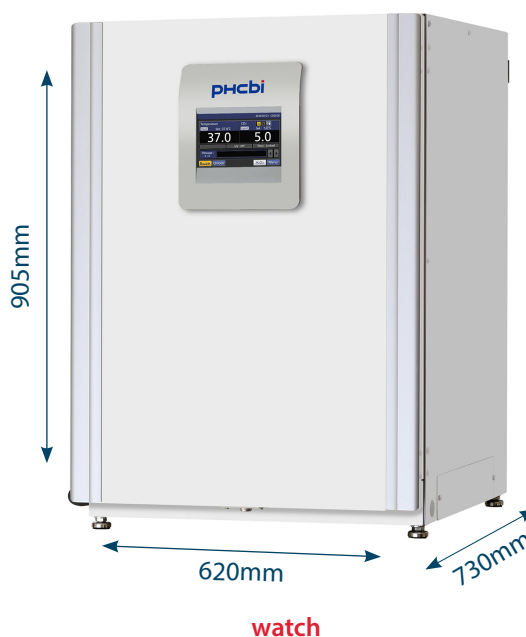

Prices and discounts may be subject to change and only available while stocks last. Not valid with any other special offer or discounts and cannot be applied in retrospect. All prices in local currency excluding GST and local delivery charges. Products advertised may be for research use only (RUO) or in vitro diagnostic use (IVD). Please refer to manufacturer's product insert for further information. For product use outside of that specified on the manufacturer's labelling and/or instructions for use, Bio-Strategy does not take responsibility.

Ref. Code:
SANY-MCO170

VALID UNTIL:
31 Mar 2024