


Energy Efficient ULT 729L Freezer

Only \$17,528

Was \$25,040

Code: SANYMDF-DU702VH-PE

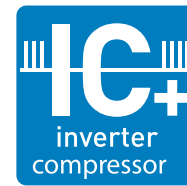

The **VIP ECO Green Gas -86°C** ultra-low temperature freezer uses natural refrigerants combined with inverter technology to provide efficient cooling and lower operating costs.



Storage capacity:

729 litres = 576 x 2" cryoboxes

- Quiet in operation – only 52dB(A) and runs on a 220-240V 50Hz 10 amp supply
- Alarms for high/low temperature, door open, power failure and filter
- Automatic vacuum relief port for quicker access to samples after door closure
- Hot gas bypass combined with fivefold door gasket assists in minimising frosting and ice build up around door seal
- Three access ports
- Three adjustable height stainless steel shelves creating four internal compartments
- Natural (HC) refrigerants, inverter compressors and integrated electronics combine to lower the freezer running costs by up to 30% over conventional models
- User friendly Large Colour Graphical Control Panel with digital display (can be used with gloved hand)
- USB port for downloading logged data
- VIP Plus insulation: provides typically 24 % more storage capacity compared to conventionally insulated freezers with the same footprint
- Castors for mobility & locating feet for working position
- Ezlatch door handle allows easy sample access and the keyed door lock provides security
- Dry remote alarm contacts
- 5-year warranty

Prices and discounts may be subject to change and only available while stocks last. Not valid with any other special offer or discounts and cannot be applied in retrospect. All prices in local currency excluding GST and local delivery charges. Products advertised may be for research use only (RUO) or in vitro diagnostic use (IVD). Please refer to manufacturer's product insert for further information. For product use outside of that specified on the manufacturer's labelling and/or instructions for use, Bio-Strategy does not take responsibility.

Ref. Code:
SANY-ECO_Green

VALID UNTIL:
31 Mar 2024