




Next Generation Incubator for Optimum Cell Culture

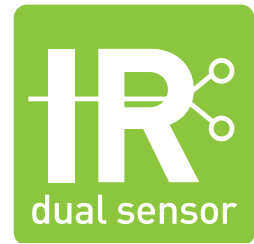

Only \$8,754

Usually \$11,672

Code: SANYMCO-230AICUVL-PE

230 litres supplied with
4 shelves & 5-year warranty

Enhance your cell growth with an intelligent CO2 incubator designed for precise temperature and CO2 control, efficient cleaning and rapid decontamination.



Stackable with
optional stacking kit

- InCusaFe™ copper alloy stainless steel interior with the benefits of copper without the tarnishing
- Temperature: Ambient +5°C to 50°C ($\pm 0.1^\circ\text{C}$)
- CO2: 0 to 20% ($\pm 0.15\%$)
- Direct heat and air jacket for environment stability
- Dual wavelength single beam IR CO2 sensor with continuous auto zero calibration
- Programmable SafeCell™ UV system (isolated from cell cultures) irradiates airflow and humidity pan to prevent contamination without affecting cell cultures in-vitro
- Humidity: 95%RH ($\pm 5\%$ RH)
- Peltier powered dew stick to draw condensation away from inner glass door
- USB port for data transfer
- Fully rounded corners and integrated shelf supports reduce contamination and simplify cleaning
- Integrated side wall tray supports – no bolts or screws that can harbour contamination

Ref. Code: SANY-MCO230

Valid Until: 31 Aug 2025



0800 34 24 66
sales.nz@bio-strategy.dksh.com
www.bio-strategy.com
labshop.dksh.co.nz



Prices and discounts may be subject to change and only available while stocks last. Not valid with any other special offer or discounts and cannot be applied in retrospect. All prices in local currency excluding GST and local delivery charges. Products advertised may be for research use only (RUO) or in vitro diagnostic use (IVD). Please refer to manufacturer's product insert for further information. For product use outside of that specified on the manufacturer's labelling and/or instructions for use, Bio-Strategy does not take responsibility.