

Separator Kit For The Disposal Of Petri Dishes

The petri dish is one of the bio waste disposables (infectious sample) which is contained cultured bacteria on the media.


“Separator Kit for Disposal Petri Dish” is the easiest solution to sterilize bio waste in disposable and separate the petri dish and media in a single autoclave cycle.

This accessory works for disposable (plastic) petri dishes only. Not suitable for glass petri dishes.



Melted petri dishes remain in the long basket after the sterilization.

The ager melts and falls down into the bucket.

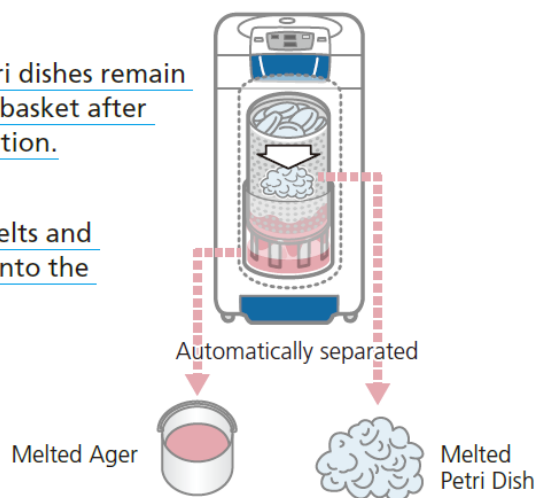

[Watch >](#)




ST-ZERO Bio waste Sterilisation Container

Bio waste disposables (infectious sample) should be sterilized before removing from the user's facility.

How do users carry these samples from bio safety cabinet to Autoclave? ST-ZERO is one solution for this issue and can be used with all models of TOMY Autoclaves.

Designed for liquid containment with flat placement and vertical mounting. Air flow will not be disturbed inside the work area of bio-safety cabinet.

- easy installation
- capacity for up to 25 x 25ml disposable pipettes
- leak proof while sterilising in vertical autoclaves

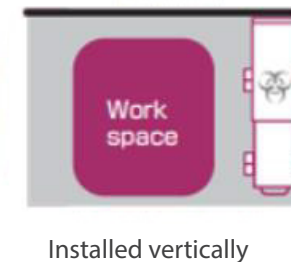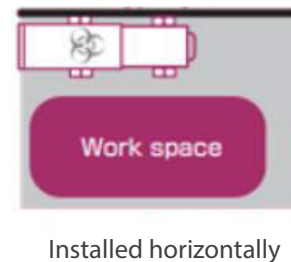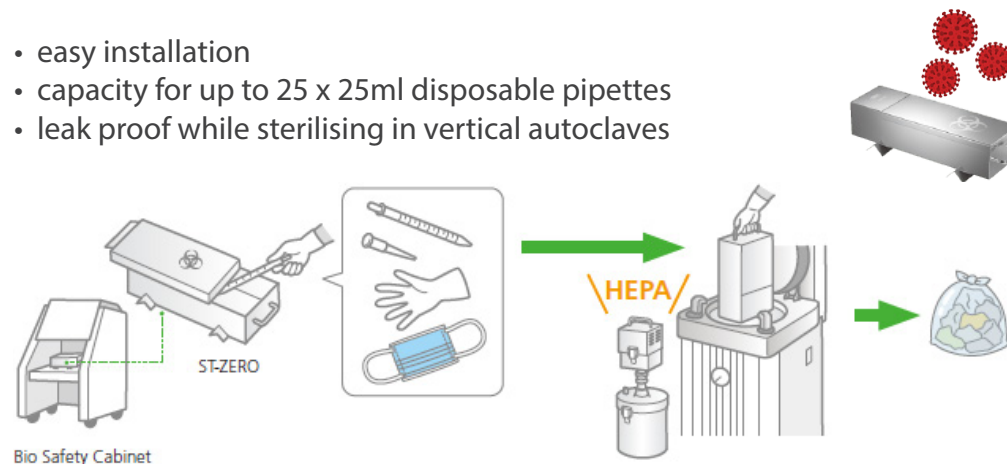

[Download manual >](#)