

PHCbi MDF-U33V -86°C Upright Freezer, 333 Litres



Valid Until:
31/12/2018

Ref. Code:
MDF-U33V

Reliable compact -86°C storage with cascade refrigeration & VIP PLUS insulation


PHCbi VIP PLUS Freezers ensures maximum storage capacity with a compact footprint.

The PHCbi MDF-U33V -86°C freezer can store 216 x 2" cryoboxes in the 333 litre capacity chamber. The vacuum panel insulation technology (VIP PLUS) and cascade refrigeration technology provide a reliable unit with increased compressor longevity and lower heat output. The unit is quiet in operation, runs off a standard 10-amp supply and the easy access airflow filter reduces maintenance times.

- Cascade refrigeration provides proven reliable cooling (-86°C in ambient of 30°C)
- **VIP Plus insulation:** provides up to 30% more storage capacity compared to conventionally insulated freezers with the same footprint
- User friendly eye-level LED display with audible and visual alarms
- Two removeable insulated inner doors with eagle-latches reduce heat loss & allow easy cleaning
- Quiet in operation: 49dB(A)
- Keyed door lock for security

In April 2018, Panasonic Healthcare Co., Ltd. changed company name to PHC Corporation. PCHbi is the new brand of PHC Corporation's biomedical division.

Speak to your local Bio-Strategy Representative for further details or call Customer Service on 0800 34 24 66 or email sales.nz@bio-strategy.com

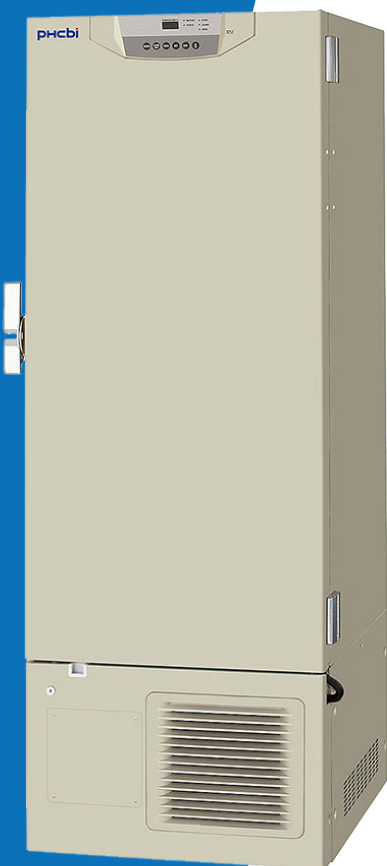

Cascade refrigeration combined with vacuum insulation panel technology

Description	Code	List Price (excl. GST)	Sale Price (excl. GST)
-86°C Upright Freezer, 333 litres	SANYMDF-U33V	\$13,000.00	\$11,675.00



SALE

BIO-STRATEGY LIMITED

Tel: 0800 34 24 66
Email: sales.nz@bio-strategy.com
www.bio-strategy.com



bio-strategy
delivering technology

Prices & discounts may be subject to change and only available while stocks last. Not valid with any other special offer or discounts and cannot be applied in retrospect. All prices exclude shipping & delivery charges. Promotion code (if stated) must be referenced on your purchase order to receive offer prices.